



Thin tube replacement – Ambio

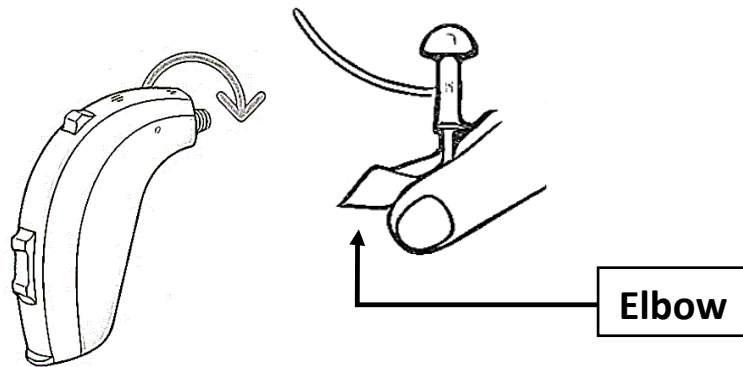

Unscrew the elbow (the plastic attachment connecting the tube and the hearing aid) anti-clockwise to remove the tube from the hearing aid. Screw the elbow of the new attachment in a clockwise direction into the hearing aid until you cannot rotate it any further. Please avoid twisting the tube while doing the replacement.

Identifying left and right

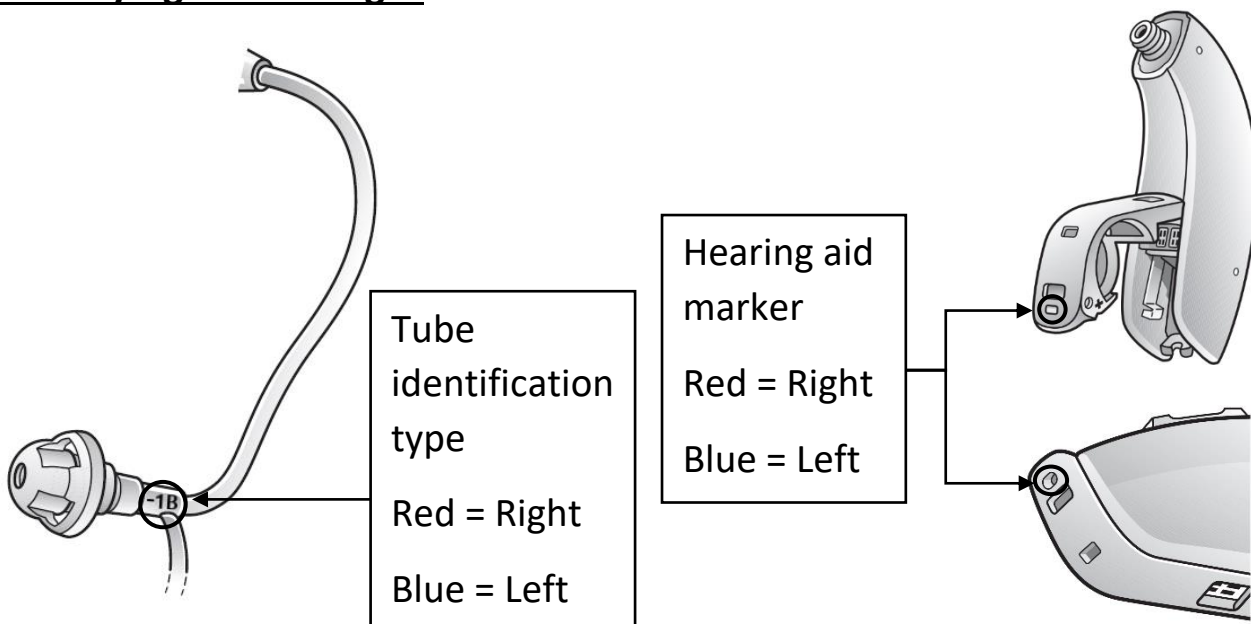

Your hearing aids are programmed to match your hearing loss for each ear individually, therefore it is important to wear the correct hearing aid on the correct side. At the bottom of the hearing aid you should see a red or blue dot (known as the hearing aid marker). The red marker indicates the right hearing aid and the blue marker indicates the left hearing aid. Likewise, at the bottom of your tube, you should see a number and a letter in red or blue. Please ensure you match the correct tube with the correct hearing aid.